# HOLOOR The Early Pioneer of Diffractive Optics

## **Vortex Lens**

Our Vortex lenses convert a Gaussian input profile into a donut-shaped energy ring.

### Features:

- High power threshold
- High efficiency
- Low back reflection
- Wavelengths from UV to IR
- Optional AR/AR coating

- Applications:
  Solar coronagraphs (astronomy)
- High-resolution microscopy
- Laser Welding
- Optical tweezers for particle trapping & manipulation
- Quantum optics

The Vortex lens is a unique optic, its structure composed entirely of spiral or helical phase steps, whose purpose is to control the phase of the transmitted beam.



The topological charge, denoted in the literature as m, refers to the number of  $2\pi$  cycles (i.e. "staircases") etched around 360° turn of diffractive surface. In Fig.1 below, the surface profiles are illustrated for vortex lenses with m=2, m=3 and m=4.

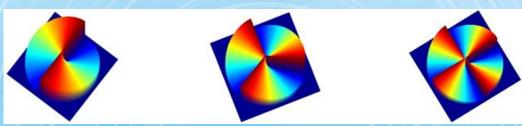


Figure 1: Surface profiles for Vortex lenses with m=2, m=3 and m=4.

One main effect of a higher topological charge is an increase in the angular momentum of the vortex beam by a factor of m. Another effect is the dimensions magnification of the ring intensity pattern, by a factor of m.



© Copyright 2018, Holo/Or Ltd. — 13B Einstein St., Ness-Ziona, 7403617, Israel Phone: +972 (8)-9409687, Fax: +972 (8)-9409606, website: <u>www.holoor.com</u>

## HOLOOR The Early Pioneer of Diffractive Optics

## **Design Considerations**

The Vortex lens requires а collimated Single Mode (TEM<sub>00</sub>) Gaussian input beam, which it converts a TEM<sub>01\*</sub> axially to symmetric mode. (\*Laguerre mode). For small divergence angles (<3°), only the working distance is affected, not the output quality.

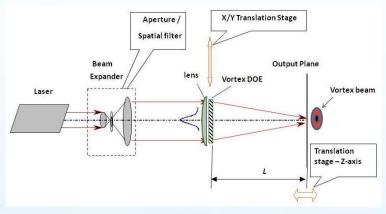


Figure 2: Typical Set-up for Vortex beam system.

The added spatial filter/aperture acts to reduce parasitic modes, and the beam expander to increase the incident beam diameter.

There are two main advantages for working with a larger input beam:

1) It reduces the sensitivity of the output to DOE alignment tolerances.

2) Enables smaller vortex spot.

The two translation stages in fig.2 above are meant to provide the user precise control of elements' locations, to reduce tolerance effects.

Spot size in the image plane of the vortex of charge m will be: (m+1) x DL, where the formula for the Diffraction-Limited spot diameter at  $1/e^2$  follows:

$$\frac{4 \times L \times \lambda}{\pi \times D} \times M^2 = DL.spotsize$$

L = Working Distance,  $\lambda$  = Wavelength, D = Input Beam Size, M<sup>2</sup> = quality of input laser beam.

### Specifications:

Materials:	Fused Silica, ZnSe, Plastics
Wavelength range:	193nm to 10.6um
DOE design:	Binary, 8-level, 16-level
Element size:	5mm to 38.1mm
Diffraction efficiency:	75%-98%
Coating (optional):	AR/AR
Topological number:	M=1, 2, 3, 4, 5, 6
Custom Design:	Ring thickness, Square donut



© Copyright 2018, Holo/Or Ltd. — 13B Einstein St., Ness-Ziona, 7403617, Israel Phone: +972 (8)-9409687, Fax: +972 (8)-9409606, website: <u>www.holoor.com</u>